



## Tenure Quik with Fluoride Universal Light-Cure, Single-Bottle Bonding Agent




REF 031145700

### CONTENTS:

- Tenure Quik with Fluoride 5 ml **CE 0044**
- Etchant Medium Viscosity 3 ml **(NM-D)**
- Syringe tip 18 G, 6 pc. **(NM-D)**
- Dappen Dish, 2 pc. **(NM-D)**
- Brush tip handle **(NM-D)**
- Brush tip, 25 pc. **(NM-D)**
- Instruction sheet **(NM-D)**
- All-Surface Bonding Chart **(NM-D)**

**Do Not Refrigerate**

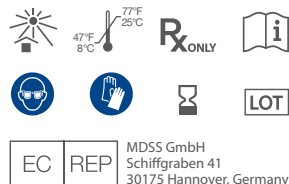

**Symbols Glossary available at:**  
[denmat.com/symbols](http://denmat.com/symbols)

Den-Mat Holdings, LLC  
1017 W. Central Avenue  
Lompoc, CA 93436 USA

©2018 Den-Mat Holdings, LLC. All rights reserved.  
751145760 R:8 11/18CM



1. Etch tooth surface with 37% phosphoric acid etchant for 15 seconds. Thoroughly rinse and gently blot or air-dry the tooth surface.



2. Completely immerse the brush tip in Tenure Quik® with Fluoride. Using a fully saturated brush tip, apply three consecutive coats to the prepared tooth. Allow Tenure Quik with Fluoride to sit on the tooth for 15 seconds.

3. Remove the solvent from the Tenure Quik with Fluoride by directing a gentle stream of dry, oil-free air from the air syringe. The force of the air should be strong enough to extract the adhesive's solvent, but not so forceful as to push the adhesive from the preparation. If the surface is not glossy after drying, reapply Tenure Quik with Fluoride and gently re-dry.



4. Light-cure Tenure Quik with Fluoride starting with:
  - Sapphire PAC lights (all models) start with 5-second exposures.
  - Flashlite LED lights (all models) start with 10-second exposures.
  - For all other curing lights (halogen, LED and other) refer to the manufacturer instructions. A minimum of 10-30 seconds is recommended.

Apply restorative material as directed by the manufacturer.

**CAUTION:** Wear protective gloves while using this product.

**CAUTION:** Wear eye protection while using this product.

**CAUTION:** Tenure® has not been studied in children, pregnant or breast-feeding women.

#### **Contraindications**

Patients allergic or hypersensitive to methacrylate monomers.

#### **INSTRUCTIONS FOR ROOT/DENTIN DESENSITIZATION**

1. Clean the sensitive area with Porcelain Polishing Paste.
2. Rinse with warm water and blot dry with a cotton pellet.
3. Apply multiple coats (4-6) of Tenure Quik with Fluoride to the surface. Air-dry with a syringe for several seconds to remove the solvent. Light-cure Tenure Quik with Fluoride as directed in Step 4.
4. Repeat Step 3.
5. For longer lasting desensitization, apply a thin layer of restorative material (for example, Virtuoso® Flowable) on top of the Tenure Quik With Fluoride.

#### **STORAGE**

**Do not** freeze.

**Do not** refrigerate.

**Do not** expose to temperatures exceeding 77° F (25° C).

**Do not** expose to direct sunlight.

Recommended temperature when used is 70° – 75° F (21° – 24° C).

Cool storage is recommended in tropical climates.

Bring to room temperature prior to use.

**SAFETY DATA SHEETS AVAILABLE AT [denmat.com](http://denmat.com)**