

## CeriSaw Posterior Finishing System for Posterior Use



REF 031336600



### CONTENTS:

- CeriSaw Posterior with 1 blade **CE**
- CeriSaw Posterior blades, 10 pc. **CE**
- All-Surface Bonding Chart **(NM-D)**
- Instruction card **(NM-D)**



Symbols Glossary available at:  
[www.denmat.com/symbols](http://www.denmat.com/symbols)

EC REP MDSS GmbH  
Schiffgraben 41  
30175 Hannover, Germany  
©2018 Den-Mat Holdings, LLC. All rights reserved.  
753366150 R:7 06/18CM

Den-Mat Holdings, LLC  
1017 W. Central Avenue  
Lompoc, CA 93436 USA

The CeriSaw is a hand-held, interproximal saw specially designed to open fused contacts during bonding, with minimal force using only one hand. The CeriSaw is available in two sizes: anterior (regular size) for anterior interproximal surfaces, and posterior (small size) for posterior interproximal surfaces. The saw is available in kits with 1 handle and 11 blades. Replacement blade packages contain 10 blades. The serrated edge of the CeriSaw blade is approximately 0.5-0.7 mm thick. The CeriSaw can be sterilized in an autoclave.

### CAUTION

The CeriSaw uses high quality, surgical blades. Extreme care should be taken when using these sharp instruments to avoid injuring the interdental papilla. By gently placing an interproximal wedge, a safety stop is provided to help prevent accidental cutting of the papilla.

Note: When using with Lumineers®, wait a minimum of 24 hours before opening contacts to ensure maximum bonding has occurred.

### DIRECTIONS

1. Completely cure the composite restorative.
2. Thoroughly remove gross excess from the interproximal regions.
3. Place the CeriSaw at the interproximal area you want to open.
4. Saw with a rocking motion and continue until the entire contact has been opened.
5. Complete normal finishing procedures.



It is important to keep the CeriSaw blade taut by rotating the handle grip sleeve. Do not over tighten the blade. Reverse the rotation of the handle grip sleeve for blade replacement.

### Adding a strip

1. Loosen the handle grip by turning handle grip counter-clockwise.
2. Wrap the CeriSaw strip onto the pin located at the front end of the handle, then back around to the pin on the rear end of the handle. The strip forms an "S" shape around the front and back handles. (Fig. 1)
3. Turn the handle grip clockwise to tighten the handgrip, which also tightens the pins and CeriSaw strip. (Fig. 2)

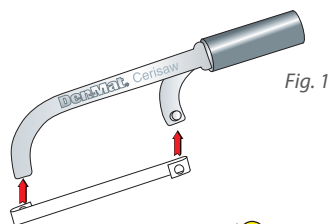


Fig. 1

### Removing a strip

1. Loosen the handle grip by turning the handle grip counter-clockwise.
2. Lift the CeriSaw strip from the pins.

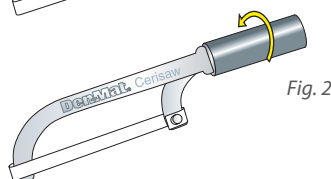


Fig. 2

## RELATED PRODUCTS AVAILABLE FROM DEN-MAT:

DESCRIPTION	PRODUCT NUMBER
CeriSaw™ Anterior, Kit ..... <i>CeriSaw Anterior with Blade, 1 ea.</i>	031336510
CeriSaw Posterior, Kit..... <i>CeriSaw Posterior with Blade, 1 ea.</i>	031336610
CeriSaw Anterior, Kit..... <i>CeriSaw Anterior with Blade, 1 ea.</i> <i>CeriSaw Anterior Blades 10, ea.</i>	031336500
CeriSaw Posterior, Kit..... <i>CeriSaw Posterior with Blade, 1 ea.</i> <i>CeriSaw Posterior Blades, 10 ea.</i>	031336600
CeriSaw Anterior, Blades, 10 ea. ....	031336520
CeriSaw Posterior, Blades, 10 ea. ....	13371

## LIMITED WARRANTY

DenMat will replace, or refund the purchase price of any of its products that are proven to be defective within 30 days of the purchase date. Replacement of defective goods, or refund of the purchase price shall be the exclusive remedy for the user. DenMat will not be liable for any economic, incidental, or consequential loss or damage that arises out of the use of, or the inability to use its products.

This limited warranty is in lieu of all other warranties, expressed or implied, and shall be void if the product is improperly stored or used. There are no implied warranties of merchantability, fitness for a particular purpose, or otherwise. Before using this product, the user shall determine whether it is suitable for the intended use and the user shall assume all risk and liability associated therewith.

## DENMAT CUSTOMER SUPPORT

U.S. and Canada.....	800-433-6628
U.S. and Canada Fax Orders and Inquiries.....	800-233-6628
International Inquiries.....	805-347-7990
International Fax Orders and Inquiries: .....	805-347-7983
Website .....	www.denmat.com

**SAFETY DATA SHEETS AVAILABLE AT [WWW.DENMAT.COM](http://WWW.DENMAT.COM)**